




CORDLESS WATERFLOSSER

User Manual

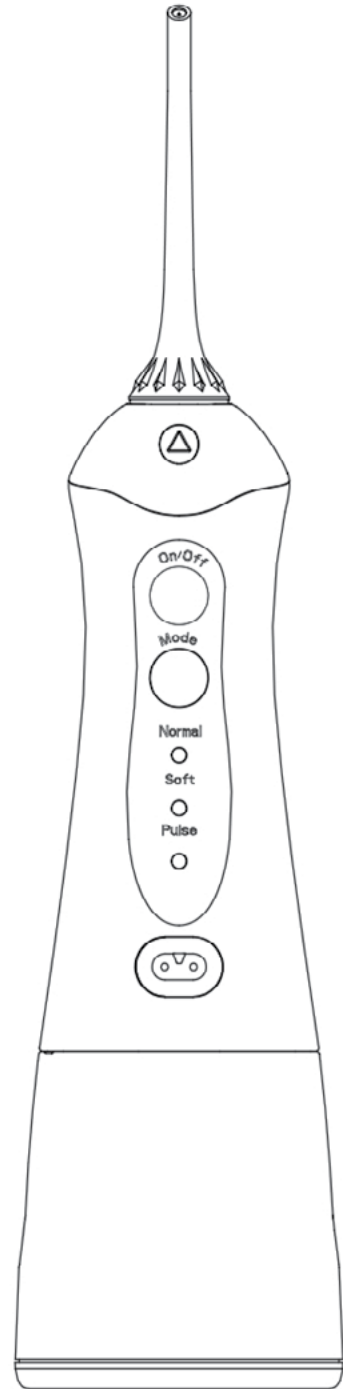

CONTENTS

| | |
|--------------------------------|----|
| SAFETY INFORMATION | 03 |
| DANGER | 03 |
| WARNING | 03 |
| IMPORTANT INFORMATION | 04 |
| PRODUCT DESCRIPTION | 05 |
| PRODUCT ACCESSORIES | 06 |
| OPERATING PROCEDURES | 07 |
| INSTRUCTIONS FOR USE | 08 |
| RECOMMENDED TECHNIQUES | 09 |
| DISASSEMBLY EXAMPLE | 10 |
| BATTERY MAINTENANCE & HANDLING | 11 |
| COMMON PROBLEMS AND SOLUTIONS | 12 |
| 360PRO THANK YOU/CONTACT | 13 |

SAFETY INFORMATION

Before using this product, please read this user manual carefully and store for future reference.

DANGER

- Do not handle charger with wet hands.
- Do not immerse in water or other liquid.
- Do not use in the bath or shower.
- Do not place or store product where it can fall or be pulled into a bath shower or sink.
- Do not reach for a product that has fallen into water. Unplug immediately.
- Check the charger cord for damage before the first use and during the life of the product.

WARNING

- Children should be supervised to ensure that they do not play with the appliance.
- The appliance is not to be used by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
- WARNING: For the purposes of recharging the battery, only use the detachable supply unit provided with this appliance.
- This symbol refers to detachable supply unit:




Model: ES050C-U050100XOS

- This appliance contains batteries that are non-replaceable.


IMPORTANT INFORMATION

- A product should never be left unattended when plugged in.
- Close supervision is necessary when the product is being used by children.
- Use this product only for its intended use as described in this manual.
- Do not use attachments not recommended by the manufacturer.
- Never operate this product if it has a damaged cord or plug. If it is not working properly or is damaged, or has been dropped into water, return the product to the 360PRO Service Centre for examination and repair.
- Keep the cord away from heated surfaces.
- Never block the air openings of the product or place it on a soft surface such as a bed or couch. Keep the air openings free of lint and hair.
- Do not use outdoors or operate where aerosol products are being used or where oxygen is being administered.
- Connect this product to a properly grounded outlet only.
- Fill reservoir with water or diluted mouth-rinse only. Never use salt or baking soda.
- In addition this user manual is also available on our website: www.my360pro.co.nz

PRODUCT DESCRIPTION

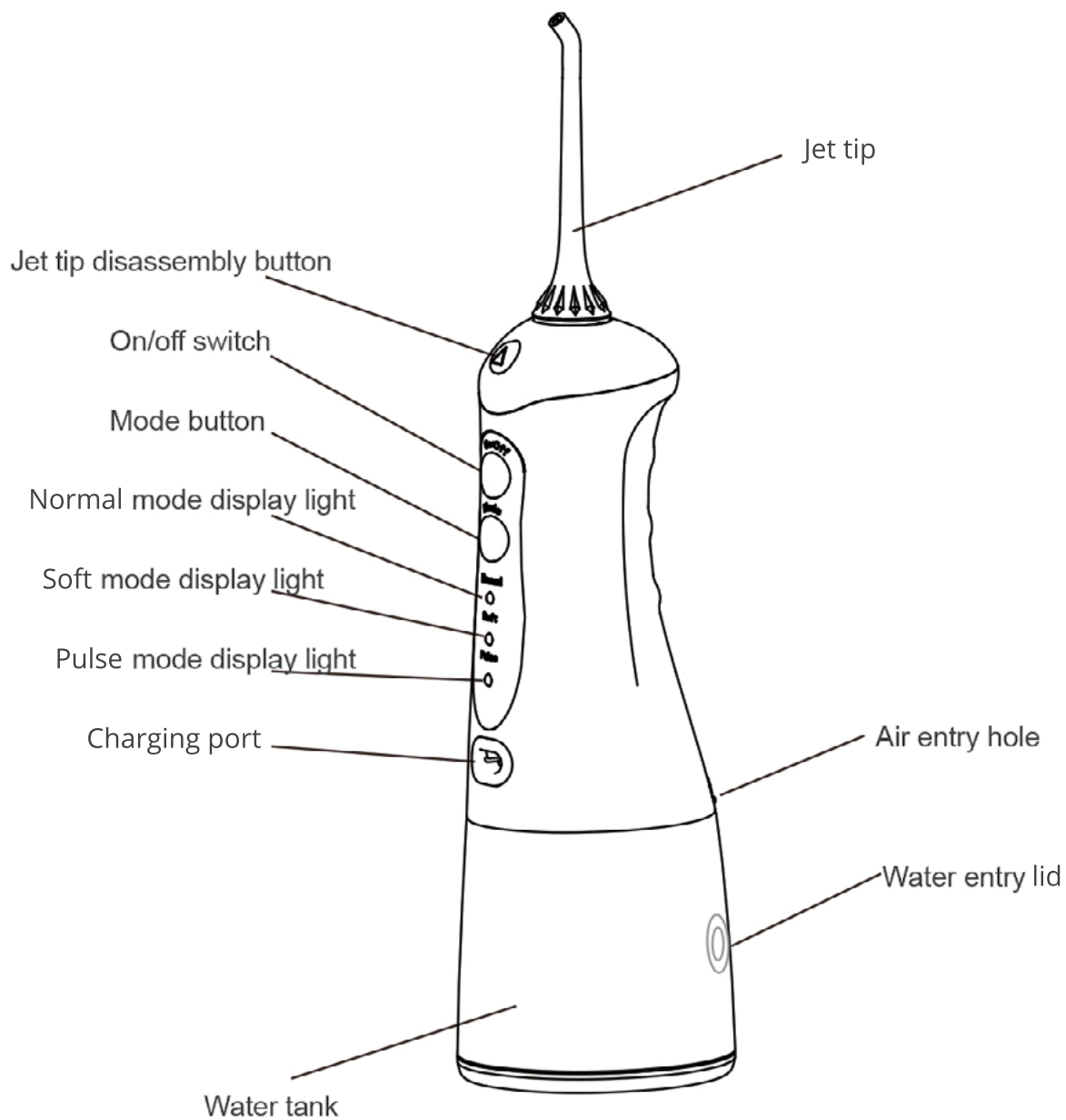

PRODUCT ACCESSORIES

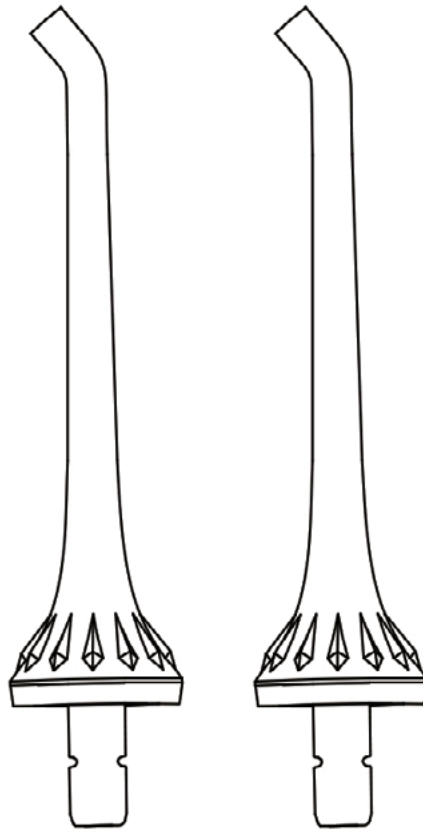

Standard jet tips x 2

The jet tips and other changable components can be purchased at www.my360pro.co.nz

It is recommended that you replace your standard jet tips every 3 months or less if they show signs of wear and tear.

OPERATING PROCEDURES




SET-UP

Extremely important!!


Prior to using for the first time, please fully charge the device (for at least 4 hours), and insert the plug into the power socket, as well as connect to a power supply. During charging, the product cannot be used.

INSTALLATION AND DISASSEMBLY OF TIPS


Insert the standard jet tip into the spray head hole until you hear a click.

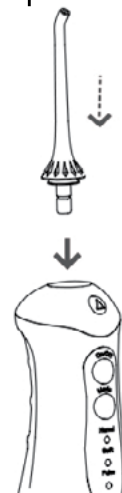

Open the water tank lid and pour water into the tank. Note: please ensure that when pouring water into water tank, the device is switched off.

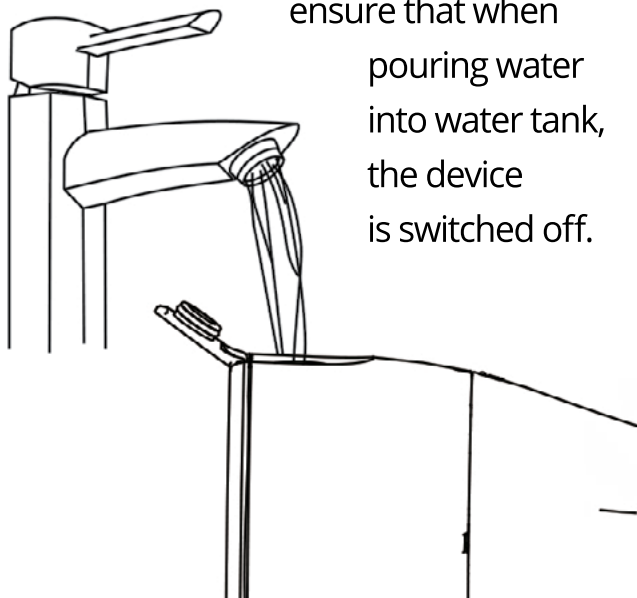

If changing or removing the tip, please press on the disassembly button on the jet tip handle before removing the jet tip.

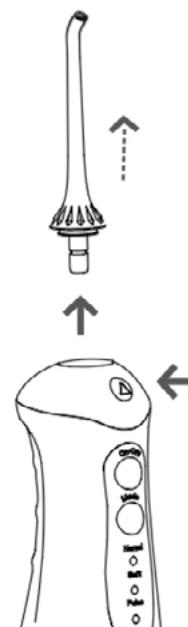


INSTRUCTIONS FOR USE

- The unit has 3 modes: Normal, soft and pulse.
- The modes may be changed when the device is switched off, and you can select the mode you like before switching the device back on.
- Irrespective of the mode, after operating for 120 seconds the product will automatically switch off.
- During recharge, the three display lights will flash, when all three lights are constantly on, this indicates that the unit is recharged.
- During usage, if all three display lights are flashing, this indicates that the products power is insufficient, and recharging is required.
- During recharge, this product cannot be used.
- Rotate the jet tip dial to adjust the direction of the water flow.




RECOMMENDED TECHNIQUES

Holding the device at a downward slant, position the jet tip in the mouth, and direct the jet tip at teeth before switching on the device.

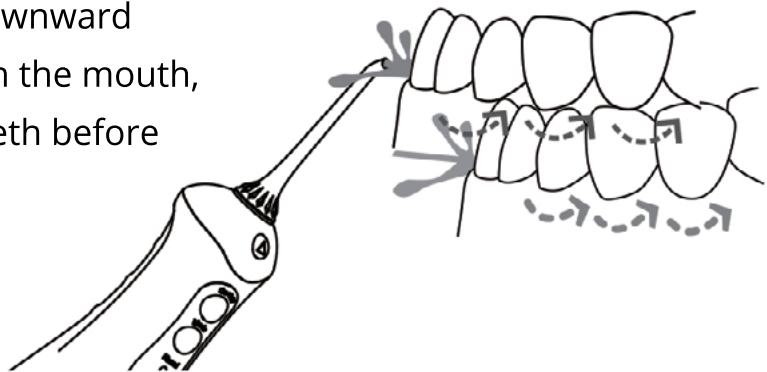


When the jet stream is at 90 degree angle to the gum line, slightly close your mouth to avoid the water spraying out. Slide the tip along the gum line as well as pausing briefly between teeth until the area around the gums and the spaces between the teeth are fully cleaned.

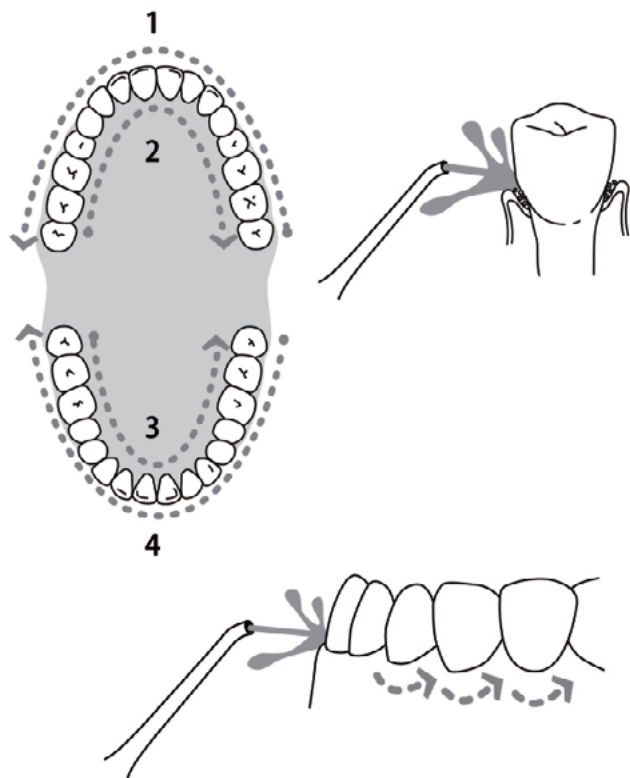


DISASSEMBLY EXAMPLE

AFTER USE

Please switch off the power, as well as emptying the water tank.

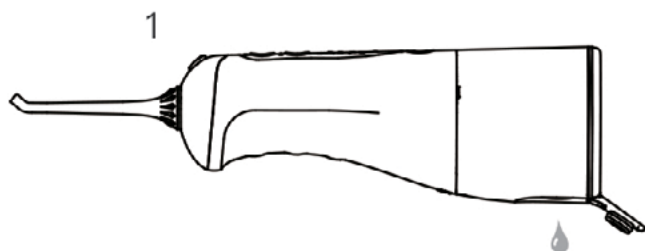

Using both hands as indicated in Diagram 2 – use your left hand to hold the product body, and your right hand to hold the water tank. At the same time rotate in an anti-clockwise direction, in order to remove the water tank. Empty.

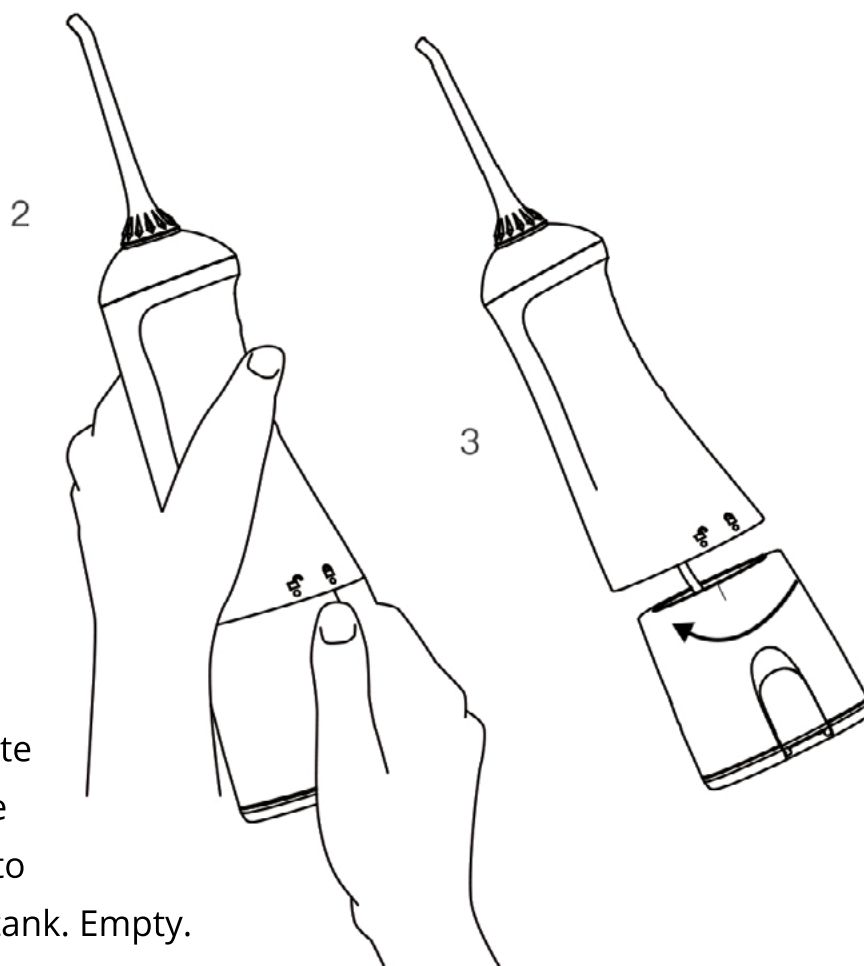

BATTERY MAINTENANCE & HANDLING

BATTERY MAINTENANCE

- Prior to initial usage the device must be fully charged, for approximately four hours.
- If to be stored for a long period of time, please fully charge the device in advance.
- In order to maximize the usage life of the battery, please promptly recharge when power is insufficient.

BATTERY HANDLING

- The device includes a lithium battery which cannot be changed. Towards the end of the battery life, please dispose of responsibly.

COMMON PROBLEMS AND SOLUTIONS

| Problem | Reason | Method of Handling |
|------------------------------|--------------------|---------------------|
| Water tank leaking water | Seal ring is loose | Reinstall correctly |
| Inadequate pressure | Low Voltage | Recharge |
| Device cannot be switched on | Low Voltage | Recharge |

If the above methods cannot resolve these problems, please consult the 360PRO Service Centre at www.my360pro.co.nz

The warranty covers defects in materials and workmanship for 2 years from the date of purchase.

Items not covered under warranty include tips, reservoir and the hose and handle.

Proof of purchase may be required.

For new tips go to www.my360pro.co.nz



THANK YOU FOR PURCHASING YOUR 360PRO CORDLESS WATERFLOSSER

Waterflossers pulse water to clean the teeth, and remove harmful bacteria below the gumline, places which cannot be cleaned by brushing the teeth alone. Manual toothbrushes are unable to effectively clean all areas of the mouth. Dental floss can only make contact with a limited area beneath the gum. The 360PRO Cordless Waterflosser can clean out harmful bacteria hidden between teeth and under the gum line. Replacement tips can be purchased at www.my360pro.co.nz



360PRO | www.my360pro.co.nz

P(+64) 9 477 0285 | PO Box 101613, NSMC, Auckland NZ 0745

Designed in New Zealand, made in China