




PLUG-IN WATERFLOSSER

User Manual

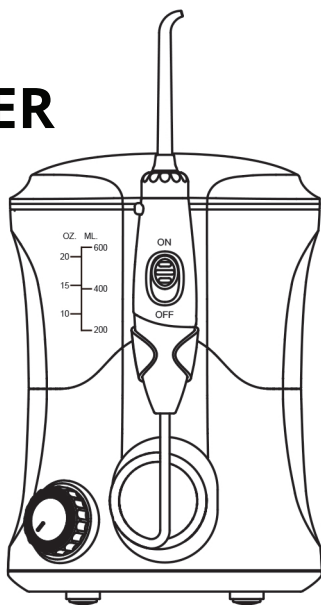

THANK YOU FOR CHOOSING YOUR 360PRO PLUG-IN WATERFLOSSER


You are on your way to better gum health and a healthier smile!
With the combination of water pressure and pulsation, the waterflosser removes harmful bacteria deep between teeth and below the gumline where brushing alone won't reach.

Traditional Clean

Traditional brushing can't easily reach hidden food debris and harmful bacteria.




String floss and interdental brushes have limited access below the gumline.

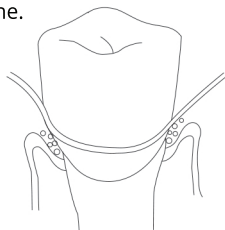

360PRO Waterflosser

Removes trapped food debris and harmful bacteria with the right combination of water pressure and pulsation.

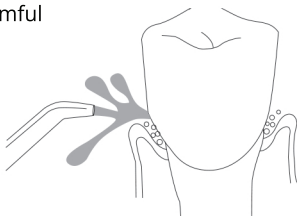

ACCESSORIES:

Main unit (including reservoir) x 1

Instruction Manual x 1

TIPS:

High-Pressure Tip x 3

Periodontal Tip x 1

Orthodontic Tip x 1

Dental Plaque Tip x 1

Tongue Cleaner x 1

IMPORTANT SAFEGUARDS

When using electrical products, basic safety precautions should always be followed, including the following:

READ ALL INSTRUCTIONS BEFORE USING

DANGER – To reduce the risk of electrocution:


- Always unplug this product immediately after use.
- Do not use while in the bath or shower.
- Do not place or store product where it can fall or be pulled into a bath, shower or sink.
- Do not place the unit or drop into water or other liquid.
- Do not place the unit or the tips in the dishwasher.
- Do not reach for a product that has fallen into water.
Unplug immediately.

WARNING – To reduce the risk of burns, electrocution, fire, or injury to persons.

- A product should never be left unattended when plugged in.
- Close supervision is necessary when the product is being used by children.
- Use this product only for its intended use as described in this manual.
- Do not use attachments not recommended by the manufacturer.
- Never operate this product if it has a damaged cord or plug. If it is not working properly or is damaged, or has been dropped into water, return the product to the 360PRO Service Centre for examination and repair.
- Keep the cord away from heated surfaces.
- Never block the air openings of the product or place it on a soft surface such as a bed or couch. Keep the air openings free of lint and hair.
- Do not use outdoors or operate where aerosol products are being used or where oxygen is being administered.
- Connect this product to a properly grounded outlet only.
- Fill reservoir with water or mouth-rinse only.
Never use salt or baking soda.

SAVE THESE INSTRUCTIONS

PRODUCT DESCRIPTION

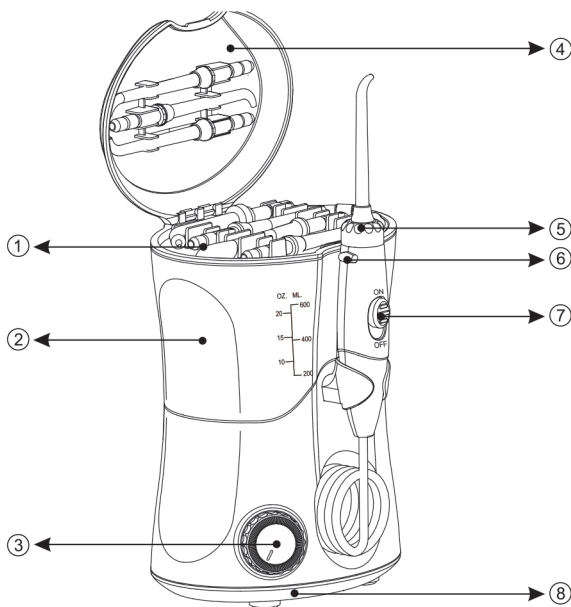

1. 7 tips

2. Reservoir

3. OFF / ON






4. Reservoir lid

5. 360° rotating tip

6. Tip removal knob

7. Pause button

8. Base suction cups

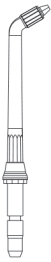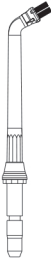				
Jet Tip	Periodontal Tip	Orthodontic Tip	Dental Plaque Tip	Tongue Cleaner


- Ordinary High Pressure Tips x 3. For every day use. Replace every 3 months.

The other 4 special tips

- Periodontal tip x 1 Replace every 3 months
- Orthodontic tip x 1 Replace every 3 months
- Dental Plaque tip x 1 Replace every 3 months
- Tongue cleaner x 1 Replace every 6 months

Tips can be removed so that multiple family members can use the same unit.

This appliance should be kept out of reach of children at all time.

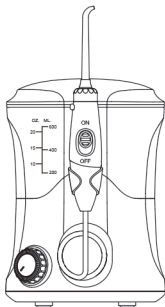


INSTRUCTIONS

Plug into Socket. If socket is controlled by a wall switch, make sure the switch is turned on.


PREPARING THE RESERVOIR

- Remove the lid & reservoir from the base.
- Fill the reservoir with lukewarm water and a small amount of mouth-rinse if required.
- Put the reservoir back on the base and press down firmly and put on the lid.

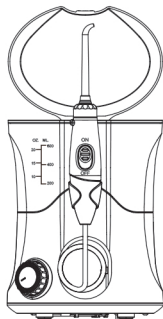



INSERTING AND REMOVING TIPS

- Insert the tip into the centre of the knob at top of the dental jet handle. The coloured ring will be flush with the end of the knob if the tip is correctly locked into place.
- To remove the tip from the handle, press the jet tip eject button and pull the tip from the handle.

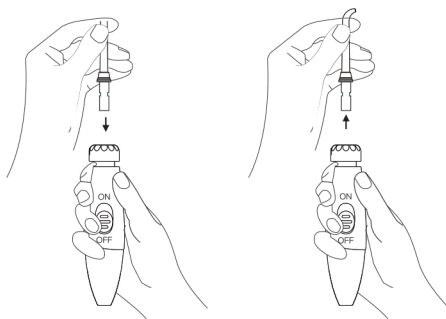




USING THE WATERFLOSSER

- When you hear click, the switch is turned on. Turn the on/off button clockwise to adjust the pressure. The machine will stop automatically after 3 minutes.

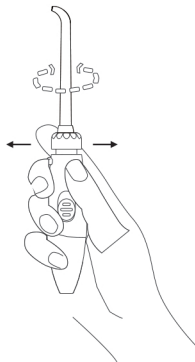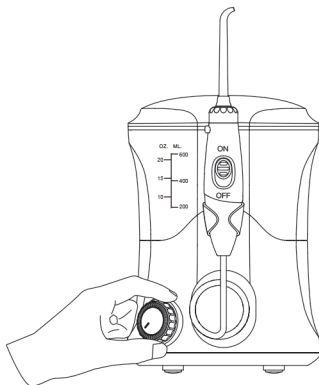


RECOMMENDED TECHNIQUE

- Direct the jet stream 90-degree angle to your gumline, slightly close your lips around the tip to avoid splashing but allow water to flow freely from your mouth into the sink.
- For best results, start in the molar area (back teeth) working towards the front teeth, glide the jet tip along gumline and pause briefly between teeth.
- Continue the process until all areas around and between the teeth have been cleaned.

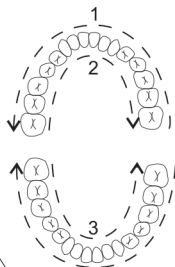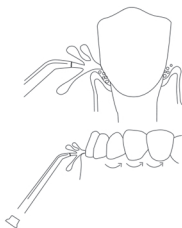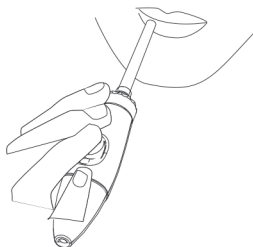


- You can temporarily stop the flow at any time by sliding the button down on the handle.
- Turn the unit off by sliding the button down and remove any liquids left in the reservoir. Unplug the unit
- The unit can also be turned off on the unit itself by turning the on/off knob anti-clockwise.

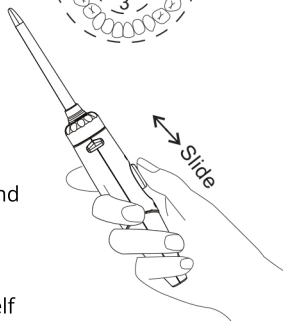

SPECIAL TIPS INSTRUCTION

PERIODONTAL TIP


- The Periodontal Tip delivers water or anti-bacterial solutions if added deep into periodontal pockets.
- To use the Periodontal Tip, lean over your sink and place the tip against a tooth at a 45-degree angle and gently place the tip under the gumline, turn the unit onto the lowest pressure setting and clean gently following the gumline, inserting tip in between teeth.

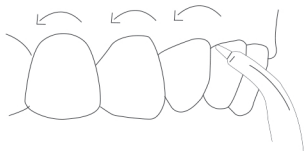

- Continue tracing along the gumline until you are finished cleaning.


TONGUE CLEANER

- Lean over the sink, place the tongue cleaner on your tongue and start the unit, use the lowest pressure setting, increase as comfortable.
- Gently place the tongue cleaner in the centre of your tongue about half way back and with light pressure pull forward.

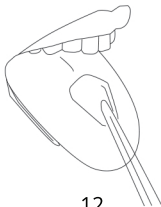


ORTHODONTIC TIP

- The special Orthodontic Tip is uniquely designed to simultaneously brush and rinse hard to reach areas around braces and other dental work.
- To use this lean low over the sink and use the brush bristles to clean around dental work. The water will flow through the brush to rinse away any debris that is removed.

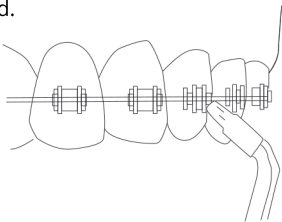

DENTAL PLAQUE TIP

- Direct the Jet Stream at a 90-degree angle to your gumline, close lips slightly to avoid splashing but allow water to flow freely from the mouth into the sink. Start in the molar (back teeth) area working toward the front teeth, place the tip close to the teeth so that the bristles are gently touching the teeth.
- Gently glide tip along gumline, pausing briefly between teeth to allow the water to flow between the teeth.
- Continue until all areas around and between the teeth both inside and outside, top and bottom have been cleaned.

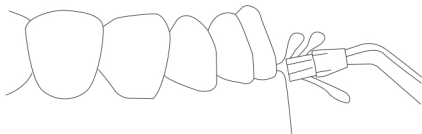

CLEANING

- After using any mouth-rinse solution in the reservoir, rinse the unit after every use.
- Partially fill the reservoir with warm water, aim the jet tip into your sink, turn on the unit and let it run until the reservoir is empty. This will help keep the tip from clogging and help your unit last longer.
- Clean the product when necessary by using a soft cloth and mild non-abrasive cleanser.
- Before cleaning, unplug from the electrical outlet.
- Regular cleaning will extend the life of your 360PRO Waterflosser.

REGULAR MAINTENANCE

- After using the 360PRO Waterflosser, unplug the power cord, and place the tank upside down on a flat surface.
- Do not disassemble or attempt to resolve any problem with the unit or the warranty will be voided.
- When the power cord is damaged, do not replace the power cord or plug yourself.
- If the product is out of the warranty period, we can provide a repair service, you will be charged the cost of materials, labour and freight if the problem can be resolved.
- The unit must only be repaired by the 360PRO service team or the warranty will be voided.

COMMON PROBLEMS AND SOLUTIONS

Problem	Cause	Solution
Leakage between the tip and handle knob	Tip not fully engaged	Remove tip and reinstall - ensure it is fully inserted
Periodontal tip - tear in soft tip	Tip is worn	Replace with new tip
Inadequate pressure	Reservoir not positioned properly	Put the reservoir on the base and push down firmly

The warranty covers defects in materials and workmanship for 2 years from the date of purchase.

Items not covered under warranty include tips, reservoir and the hose and handle.

Proof of purchase may be required.

For new tips go to www.my360pro.co.nz

Thank you again for choosing the 360PRO Waterflosser.

Thank you for purchasing your 360PRO Plug-in Waterflosser.

Waterflossers pulse water to clean the teeth, and remove harmful bacteria below the gumline and places which cannot be cleaned by brushing the teeth alone. Manual toothbrushes are unable to effectively clean all areas of the mouth. Dental floss can only make contact with a limited area beneath the gum. The 360PRO Waterflosser can clean out harmful bacteria hidden between teeth and under the gum line.

Replacement tips can be purchased at www.my360pro.co.nz



360PRO

PO Box 101 613, NSMC, Auckland, NZ 0745

P: +64 9 477 0285 | www.my360pro.co.nz

Designed in New Zealand, made in China.

